


Instructions for Oral Fluid Sample Collections with Quantisal™

Patient Self-Collection Instructions for physician-ordered tests*: For observed collection in instances when in-person collections are not available

Self-Collection Prep—Review the testing procedures with your healthcare provider, select a clean and dry place to conduct the provider-observed collection


- Check device package for integrity and appropriate expiration date. Empty your mouth of gum, food, tobacco prior to test. If mouth is not empty, rinse with water and wait 10 minutes.
- Move tongue side to side to accumulate saliva prior to collection.
- Peel open the Quantisal™ collection device package just far enough for easy removal of the collector handle. Keep the tip of the device pointed down. (Note: **Do not** use the collection pad if the package has been previously opened, if the collection pad is wet, or if indicator dye is visible prior to use.)
- Use the handle to remove the device from the packet. (Note: **Do not** touch the collection pad with fingers before or after specimen collection.)

Step 1 Position the collector under tongue and close mouth. Keep head tilted downward to allow gravity to help with saliva collection. **Do not chew on pad, talk, or remove collector from mouth until indicator turns BLUE or until 10 minutes have passed.**

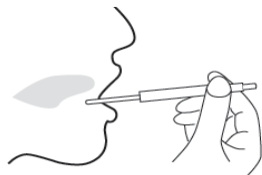


(Figure 1)

Step 2 Start a 10-minute timer. **The collection process is complete when the blue color appears in the indicator of the handle (typically 2-10 minutes), or when collection time reaches 10 minutes, whichever occurs first.**

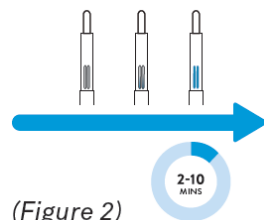


(Figure 2)

Step 3 When the collection is completed, remove the oral fluid collection transport tube from the package. Hold transport tube in an upright position and uncap by pushing up with thumb(s). **Do not** stand tube on table. **Do not** spill or empty the liquid from the tube.

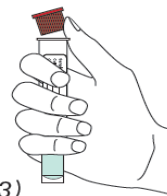


(Figure 3)

Step 4 Remove the Quantisal oral fluid collector from your mouth and guide the collector into the uncapped transport tube and replace cap. (Note: Do not place collector back in mouth after it has been placed in the oral fluid collection transport tube.)

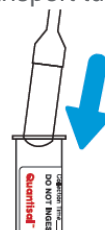

(Figure 4)

Step 5 Once the collector is inside the transport tube replace the cap. **SNAP** cap firmly in place.

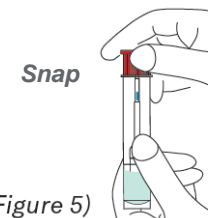

(Figure 5)

Step 6 Place the specimen vial into the specimen bag and complete the procedures to prepare and ship the specimen to the testing laboratory.

For questions about the Quantisal™ collection package, please contact your healthcare provider.

For additional training information, please visit:
www.immunalysis.com/products/oral-fluid/quantisal/collector-training/

*Intended Use: The Quantisal™ collection device is being supplied as part of a test ordered by your healthcare provider. This Quantisal™ Oral Fluid Collection Device is intended for the collection, preservation, and transport of oral fluid specimens for analytical testing of drugs or drug metabolites. This device is for use only under observed collection. These instructions are for provider-observed collections in instances when in-person collections are not available.

—Please see reverse side for packing/shipping instructions—


Instructions for Oral Fluid Sample Collections with Quantisal™

Packing & Shipping Instructions—collect specimen and ship on the same day

Step A Ensure sample collection tube is properly sealed.

Figure A

Snap

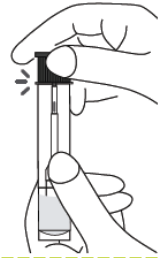

Step B Enter collection date information onto the provided requisition form.

Figure B


Collection date



Step C Peel label (detachable lower section of requisition form) and affix to sample collection tube, lengthwise down the center of collection tube.

Figure C


Peel label (lower section of form)

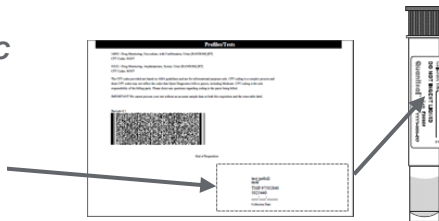

Step D Place the sample collection tube inside the provided specimen bag, in the main section of the bag.

Figure D

Tube here

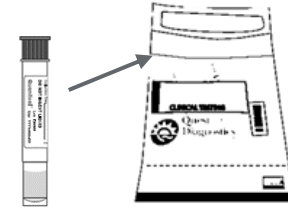

Step E Fold completed requisition form in half and place inside the specimen bag, in the smaller, front pouch section of the bag.

Figure E

Fold here

Form here



Step F Firmly close the adhesive on the specimen bag and insert into the provided shipping box.

Figure F

Insert specimen bag



Step G Close shipping box with provided seal and affix provided “Return” shipping label to sealed shipping box (over existing label).

Figure G

Affix return shipping label (over existing)



Step H Place labeled shipping box into a convenient FedEx drop box or drop off at a convenient FedEx collection site.

Figure H



For questions about the Quantisal™ collection package, please contact your healthcare provider.

For additional training information, please visit:
www.immunalysis.com/products/oral-fluid/quantisal/collector-training/

Quantisal™ is a trademark of Immunalysis Corporation or their respective owners

QuestDiagnostics.com

Quest, Quest Diagnostics, any associated logos, and all associated Quest Diagnostics registered, or unregistered trademarks are the property of Quest Diagnostics. All third-party marks — © and ™ — are the property of their respective owners. © 2022 Quest Diagnostics Incorporated. All rights reserved. PP9403 11/2022

—Please see reverse side for self-collection instructions—

